

IV ADMINISTRATION USING A BREAK & SHAKE BAG AND A FLOW RATE CONTROLLER

See instructional video by clicking link below.

<https://www.landmarkinfusion.com/flow-controlled-administration-video>

GETTING STARTED

- Use a clean work area.
- Clean your hands thoroughly.

GATHER YOUR SUPPLIES

- IV Medication bag
- IV tubing with Flow-Rate Controller
- IV pole
- Flush syringes: Saline (white cap), Heparin (yellow or blue cap)
- Alcohol Wipes
- Trash container
- SASH mat guide

INSPECT MEDICATION BAG

- Check label on medication bag.
- Verify correct name of patient, correct name of medication
- Review discard date of medication.
- Medication bag should be free of leaks or discoloration. Call pharmacy for problems with medication bag or incorrect labels.



PREPARE YOUR IV MEDICATION

- Bend up then down to break the seal of the medication vial (feels like a you are breaking a glow stick).
- Squeeze the bag to force fluid into vial.
- Gently shake vial until drug is completely dissolved.
- Hold the medication bag upside down and release squeeze to drain the mixed drug from the vial back into the bag.
- You may repeat these steps until the vial is completely empty of the drug and solution (DO NOT remove the drug vial).

**See Figure 1 for more information.*

***Figure 2 gives instructions for preparing your IV medication with an alternative medication bag that may be used.*



Irving: 8200 Tristar Drive, Suite 120
Irving, Texas 75063

Fort Worth: 3515 NW Jim Wright Freeway
Fort Worth, Texas 76106

Pharmacy Phone: (866) 388-3883
Pharmacy Fax: (888) 249-3132

www.landmarkinfusion.com

PRIME YOUR IV TUBING

- Remove IV tubing from its package. Close the clamp on the tubing.
- Place the medication bag on a clean flat surface.
- Remove the tab from the medication bag.
- Remove the cover from the spike of the tubing (DO NOT touch the end).
- Spike the bag by pushing and twisting the end of the tubing into the bag. Make sure the spike reaches the medication in the bag and be careful not to puncture the bag.
- Hang the IV medication bag on the IV pole. Gently squeeze the drip chamber to fill it about ½ full of fluid (DO NOT overfill).
- Turn the IV tubing dial to “OPEN”. Open the clamp to prime the tubing (filling the tubing with medication). When the medication reaches the end of the tubing, close the clamp.
- Adjust the tubing dial to set the correct rate for your IV infusion (your rate will be on medication label).

**INFUSE MEDICATION VIA SASH PROCEDURE AS DETAILED BELOW:****1. Flush your IV line with saline (S)**

- If applicable, release the clamp on your IV catheter and/or extension.
- Remove alcohol cap from the end of your IV catheter.
- Scrub the needleless connector on your IV catheter with an alcohol wipe for at least 15 seconds and allow to air dry.
- Remove the protective cap from the saline flush syringe. DO NOT touch the syringe tip after removing the cap. Push out air and few drops of excess solution by holding syringe with the tip facing the ceiling and pushing the plunger up.
- Attach the saline flush syringe to the needleless connector on your IV catheter using a push and clockwise twisting motion until secured together.
- Flush your IV catheter using a “push/pause” method.
- Remove and discard the saline flush syringe in a trash container.
- IF resistance is met when trying to flush, verify clamp is open. IF resistance is still met, call your nurse.



Irving: 8200 Tristar Drive, Suite 120
Irving, Texas 75063

Fort Worth: 3515 NW Jim Wright Freeway
Fort Worth, Texas 76106

Pharmacy Phone: (866) 388-3883
Pharmacy Fax: (888) 249-3132

www.landmarkinfusion.com

2. Administer Medication (A)

- Scrub the needleless connector on your IV catheter with an alcohol wipe for at least 15 seconds and allow to air dry.
- Remove the protective cap from the IV tubing.
- Attach the IV tubing firmly to the needleless connector on your IV catheter using a push and clockwise twisting motion until secured together.
- Make sure your infusion dial is set to the correct rate, then open the clamp on IV tubing.
- Check your infusion consistently to make sure it is running correctly.
- Once the prescribed dose has infused, close the clamp on the IV tubing and disconnect from your IV catheter. Discard all equipment in a trash container.



3. Flush your IV line with saline (S)

- Follow step 1 above starting at the third bullet point.

4. Flush your line with heparin (H)

- Follow step 1 above starting at the third bullet point, using heparin rather than saline.
- Close the clamp on your IV catheter and/or IV extension set, if indicated.

****IF YOU HAVE TWO PORTS ON YOUR IV CATHETER...***

- *Make sure you flush each port at least once daily. Complete (SH) portion of (SASH) procedure on unused IV catheter port by flushing with saline and then heparin as indicated above.*









Irving: 8200 Tristar Drive, Suite 120
Irving, Texas 75063

Fort Worth: 3515 NW Jim Wright Freeway
Fort Worth, Texas 76106

Pharmacy Phone: (866) 388-3883
Pharmacy Fax: (888) 249-3132

www.landmarkinfusion.com

MINI-BAG PLUS system directions:
Only for Single Dose Powdered Drug Vials with 20 mm Closures
Use Aseptic Technique

| | | |
|----------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <p>1</p>  | <p>2</p>  | <p>Assembly</p> <ol style="list-style-type: none"> Remove vial cover Disinfect stopper Peel off foil cover Inspect adaptor for moisture. Discard if found Place vial upright Hold firmly Push adaptor down until vial snaps in place DO NOT TWIST Pull vial to ensure fully seated |
| <p>3</p>  | <p>4</p>  | <p>Reconstitution</p> <ol style="list-style-type: none"> Squeeze bag and check vial Use only if vial fully seated and dry Bend up then down to break seal Hold bag with vial down Squeeze solution into vial until half full Shake to suspend drug in solution Hold bag with vial upside down Squeeze bag to force air into vial Release to drain suspended drug from vial Repeat steps 5 and 6 until vial is empty of drug and solution is thoroughly mixed. Ensure drug is completely dissolved. Do Not Remove Drug Vial. |
| <p>5</p>  | <p>6</p>  | <ol style="list-style-type: none"> Remove port protector. Attach administration set per its directions. Hang container on I.V. pole and prime set per directions. Ensure that vial is empty of drug and solution. Repeat step 6 if drug and solution remain in vial. Warning: Do not use in series connections Administer medication per directions. Use within specified time for drug stability. Refer to drug package insert. |

Rx only.

To ensure safe and proper use of this device, refer to the complete instructions for use.

To Order Baxter Products, call 1-888-229-0001



Options at work.

MEDICATION DELIVERY

Baxter, Mini-Bag Plus and Committed to a Safer Healthcare Environment are trademarks of Baxter International Inc.
 Baxter Healthcare Corporation, Medication Delivery, Route 120 and Wilson Road, Round Lake, IL 60073 www.baxter.com 111094 5K 04/09



FIGURE 1
 *see steps 4-6

ALTERNATE MEDICATION BAG

Irving: 8200 Tristar Drive, Suite 120
 Irving, Texas 75063

Fort Worth: 3515 NW Jim Wright Freeway
 Fort Worth, Texas 76106


Pharmacy Phone: (866) 388-3883
 Pharmacy Fax: (888) 249-3132

www.landmarkinfusion.com




ALTERNATE MEDICATION BAG INSTRUCTIONS


Activation and Administration Set Instructions

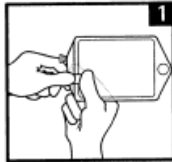


1 At time of use, check to ensure addEASE is firmly latched to container. If not firmly latched, discard assembly. Prior to reconstitution, inspect contents of drug vial for presence of liquid. If liquid is present in powdered drug, discard assembly.

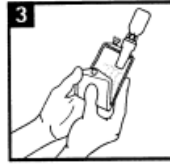


3 With vial positioned above container, as shown, fold container a third from top and squeeze firmly to transfer solution from vial to container. If necessary, repeat Steps 2 and 3.
CAUTION: Inspect solution and discard assembled system if cloudiness, haze, or particulate matter is found.
Consult medication vial Instructions for Use for additional admixing instructions.

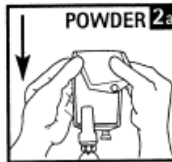

To activate vial containing *powder drug*, position vial below container **2a**. For vial containing *liquid drug*, position vial above container to activate **2b**. Fold container a third from top and squeeze firmly to force plug from spike into vial. To transfer diluent to vial, fold container a third from top and squeeze firmly. Repeat until vial is approximately half full. Do not overfill vial with liquid. Thoroughly mix to reconstitute, following drug manufacturer's Package Insert.

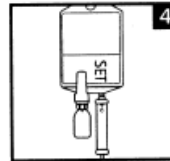

4 After reconstitution, attach the administration set to container set port. Refer to IV set Instructions for Use. Hang container. To prevent fluid from re-entering vial, avoid squeezing or flexing container.
CAUTION: Do not disconnect medication vial or PAB container from addEASE. After use, dispose of assembled system per institutional protocol.

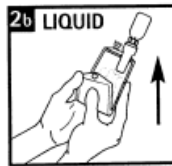

Fluid path is:

STERILE R
NONPYROGENIC

Not made with natural rubber latex or DEHP.

For single use only.
Do not resterilize.

Rx only

Compatible with most 20 mm medication vials.

Pictorials are for reference only.

B. Braun Medical Inc.
Bethlehem, PA 18018-3524 USA
1-800-227-2862



FIGURE 2

Irving: 8200 Tristar Drive, Suite 120
Irving, Texas 75063
Fort Worth: 3515 NW Jim Wright Freeway
Fort Worth, Texas 76106

Pharmacy Phone: (866) 388-3883
Pharmacy Fax: (888) 249-3132

www.landmarkinfusion.com